
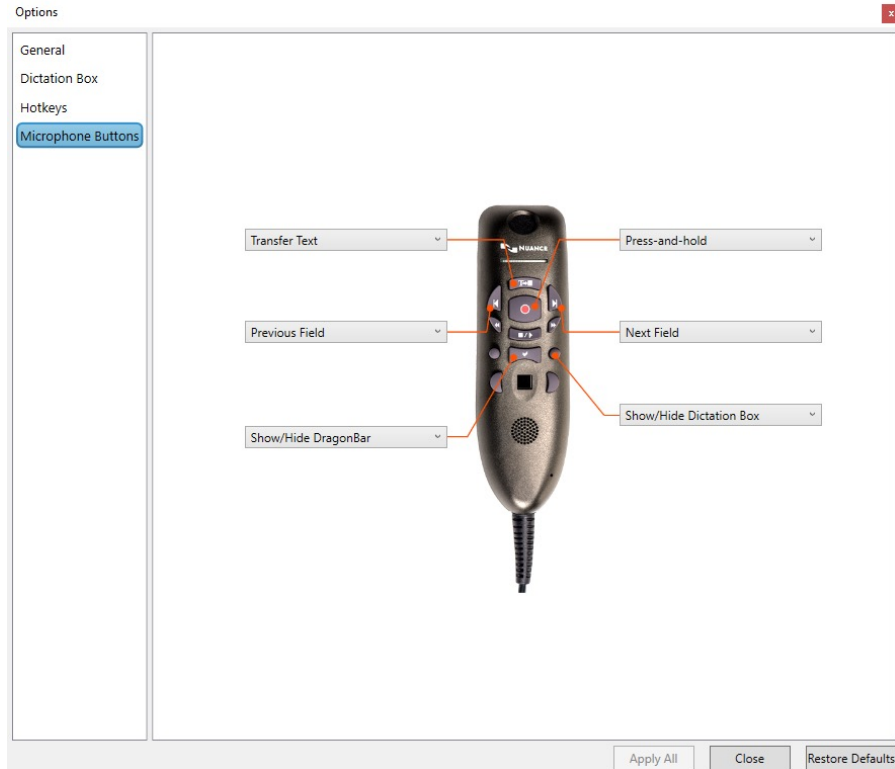


Programming the PowerMic Buttons in Dragon Medical One

To program the buttons on the PowerMic in Dragon Medical One, follow these simple steps.

1. Click the blue Dragon icon on the DMO widget, then select options.

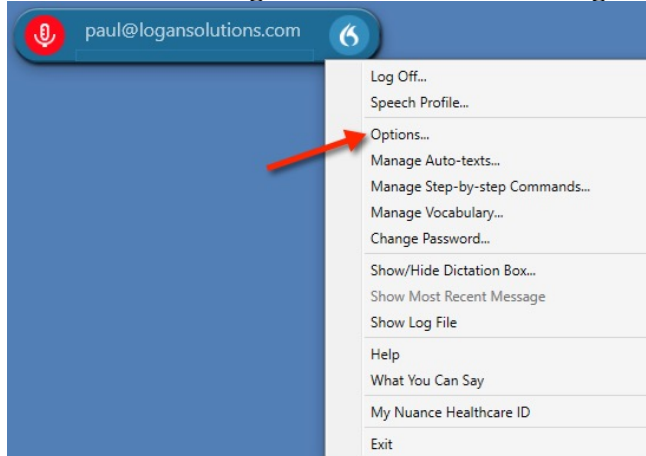

- Now, just click the button you want to change, then select the function you want for that button.

Function	Description
Clear dictation box	Removes all text from the dictation box
Recall text	Removes text transferred from the dictation box
Transfer text	Moves text from the dictation box to the document
Anchor/release speech focus	Toggles between anchoring and unanchoring the dictation to the place you leave the cursor
Show/hide DragonBar	Toggles between showing and hiding the DragonBar (e.g. widget)
Show/hide dictation box	Toggles between showing and hiding the dictation box
Next field	Moves the cursor to the next []
Previous field	Moves the cursor to the previous []
Field complete	Removes the brackets and leaves whatever text is between them
Accept defaults	Removes ALL brackets from the entire document and leaves whatever text is between them
Press hotkey	Emulates pushing a Dragon hotkey on the keyboard
Step-by-step command...	Invokes a step-by-step command you have built
Microphone on/off	Toggles between the microphone being on and off each time the button is pushed
Press and hold	Microphone is on for the entire time the button is depressed

- Click the Apply All button when you're done.